

PHILIPS

DreamStation Go

Portable PAP therapy system

Innovation for people on the Go



Convenience
without compromise



The **DreamStation Go** is a compact and highly portable PAP system designed for patients who want to maintain their therapy at home and on the road.



Powerful, portable performance

At half the size of our System One and DreamStation platforms, DreamStation Go offers the same therapy and comfort features used by more than 5 million CPAP users.



Mask compatibility

DreamStation Go gives your patients the freedom to choose the mask they want.¹



Optional battery

DreamStation Go's overnight battery can help your patients stay committed to their therapy when they are off the grid. Set at 10cm of pressure, DreamStation Go displayed an average run time of 13 hours² when operating on its overnight battery.



12mm micro-flexible tubing

DreamStation Go features the smallest PAP tubing we've ever offered, providing enhanced flexibility, portability, and connectivity. In a recent head-to-head trial, patients found the DreamStation 12mm micro-flexible tubing to offer better flexibility, more comfort, and greater ease of movement in comparison to standard tubing.³




Simple setup

With our intuitive color touch screen navigation, a mobile device app isn't required to set up or adjust your patients' therapy.



Compliance reporting made easy

Like our other PAP devices, DreamStation Go is compatible with our patient data management system, which seamlessly combines usage and compliance data from multiple Philips Respironics CPAP devices into one simple report.

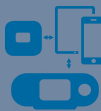

DreamMapper cross-device connectivity

Users of DreamMapper, our patient self-management tool, can sync data from multiple DreamMapper-enabled devices⁴ to one account. They also have the ability to track and monitor their combined progress across both home and portable PAP devices via one convenient dashboard on their smartphone, tablet, or computer.

¹Masks with standard 22mm outer diameter connection.

²Battery run times may vary due to higher CPAP pressure settings, excessive mask leakage, altitude, temperature, charging devices with the USB port, and other variables.

³Head-to-head trial of DreamStation 15mm tubing vs. DreamStation 12mm micro-flexible tubing, 2016.

⁴To see which therapy devices are compatible with DreamMapper, visit www.dreammapper.com/compatible.

A woman with blonde hair, wearing a blue and white striped sleeveless dress, is focused on packing a checkered bag with tan leather straps. She is standing at a table where a wide-brimmed straw hat is placed. The background is a blurred indoor setting with other people, suggesting a busy public space like an airport or train station.

The confident, **convenient choice**

DreamStation Go provides the portability and packability travelers want and the effective therapy they need.



65%
smaller*

The DreamStation Go is 65% smaller than ResMed's AirSense 10 and AirStart 10.

58%
smaller

The DreamStation Go 12mm tubing is 58% smaller than our 22mm performance tubing.

48%
lighter*

The DreamStation Go is 48% lighter than ResMed's AirSense 10 and AirStart 10.


33%
fewer components**

The DreamStation Go requires 33% fewer components to pack than its closest competitors.

*Internal measurement of ResMed's AirSense 10 and AirStart 10 with side cover installed. Includes power supply. Data on file, 2016.

**When compared to HDM Z1 and Somnetics Transcend devices, 2016.

Powerfully convenient features

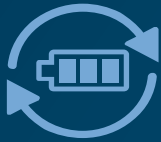

Optional battery provides approximately 13 hours of therapy!¹



Intuitive color touchscreen offers quick, simple setup and easy navigation.



Built-in power supply reduces the number of parts managed by the patient.



Optional microSD card stores compliance and waveform data.²



DreamStation Go is the only travel PAP with a USB charging port,³ allowing for easy charging of devices.



¹PAP run-time test conditions: CPAP mode, pressure 10cm H₂O, 12mm tubing, 37LPM of leak, 73.4 F room at an elevation of 1650 ft.

²Sold separately.

³When compared against HDM Z1, Somnetics Transcend, and ResMed's AirMini devices, 2017.

⁴RTCA/DO-160G section 21, category M.



Bluetooth is standard with every device.



Unit complies with FAA requirements for in-flight use.⁴



Device is compatible with DreamMapper, our patient self-management tool.

12mm tubing

Micro-flexible tubing is our smallest and lightest ever.



Click-connections

Click-connections on each end of the micro-flexible tubing make connecting and disconnecting the device and mask simple.

Create a world of opportunities



Overnight battery

Compact and lightweight, the battery features a recharge time of just 4-5 hours.



Small travel kit

Kit includes an easy-to-pack protective case for the DreamStation Go and a washable bag for the mask and tube.




Medium travel kit

Kit includes a protective case that accommodates the DreamStation Go and overnight battery, as well as a washable bag for the mask and tube.



Air filters

Choice of a reusable pollen or disposable 30-day fine filter.

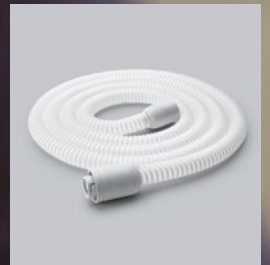

12mm micro-flexible tubing

This tubing offers more comfort and greater ease of movement than that of standard tubing.*

*Head-to-head trial of DreamStation 15mm tubing vs. DreamStation 12mm micro-flexible tubing, 2016.



DreamStation Go Auto specifications

Pressure range	4 to 20 cm H ₂ O
Pressure mode(s)	CPAP, APAP
Flex pressure relief	1 to 3
Flex modes	C-Flex, C-Flex+, A-Flex
Ramp time	0 to 45 min (5-min increments)
Starting ramp pressure	Patient adjustable, Smart Ramp
Dimensions	15 cm L x 15 cm W x 5.8 cm H
Weight	844 g (1.86 lb)
Filters	Reusable pollen, combination fine
Device setup	Touchscreen color display
Data storage capacity	On-board: > 1 year MicroSD card (optional): lifetime
Compliance meter	Breathing detection
Altitude compensation	Automatic
Electrical requirements	100-240 VAC, 50/60 Hz
Warranty	Two years

Ordering information

Flow generators	Part number
DreamStation Go CPAP	DSG400S11
DreamStation Go Auto CPAP	DSG500S11
Power	Part number
DreamStation Go overnight battery	1133281
Travel kits	Part number
DreamStation Go small travel kit*	1133275
DreamStation Go medium travel kit**	1133276
Consumables/Replacement	Part number
DreamStation micro-flexible 12mm tubing	PR12
Power cord 6ft white	1132551
Power cord 10ft white	1132552
Beauty panel - front	1133741
Beauty panel - rear	1133742
Reusable pollen filter 2-pk	1133743
Disposable 30-day fine filter 6-pk	1133744

DreamStation Go CPAP specifications

Pressure range	4 to 20 cm H ₂ O
Pressure mode(s)	CPAP
Flex pressure relief	1 to 3
Flex modes	C-Flex, C-Flex+
Ramp time	0 to 45 min (5-min increments)
Starting ramp pressure	Patient adjustable, Smart Ramp
Dimensions	15 cm L x 15 cm W x 5.8 cm H
Weight	844 g (1.86 lb)
Filters	Reusable pollen, combination fine
Device setup	Touchscreen color display
Data storage capacity	On-board: > 1 year MicroSD card (optional): lifetime
Compliance meter	Breathing detection
Altitude compensation	Automatic
Electrical requirements	100-240 VAC, 50/60 Hz
Warranty	Two years



*Small travel kit includes an easy-to-pack protective case for the DreamStation Go and a separate washable mask and tube bag.

**Medium travel kit includes a protective case that accommodates the DreamStation Go and overnight battery. Also includes a separate washable mask and tube bag.



Learn more

For more information about DreamStation Go,
visit www.philips.com/dreamstationgo.



© 2018 Koninklijke Philips N.V. All rights reserved.
Specifications are subject to change without notice.

www.philips.com/respironics

Caution: US federal law restricts these devices to
sale by or on the order of a physician.

RRDPCH EB 3/7/18 PN 1133355 MCI 4109177
1010 Murry Ridge Lane, Murrysville, PA 15668
800-345-6443 • 724-387-4000